



Sede Legale ed Amministrativa / Headquarters CEFLA s.c. Via Selice Provinciale, 23/a 40026 Imola (B0) - Italy **t.** +39 0542 653111 **f.** +39 0542 653344 sternweber@sternweber.com sternweber.com

Stabilimento / Plant

Via Bicocca, 14/C 40026 Imola (BO) - Italy t. +39 0542 653441 f. +39 0542 653601

VERSATILITY AT ITS VERY BEST

STERNWEBER.COM



S220GB191S00

04/2019











OPTIMAL FLEXIBILITY



12

Arms, modules and instruments that ensure working flexibility and comfort

Everything on the Continental model has been designed to provide the dentist with an extremely flexible system. From the compactness and lightness of the dentist's module to the swing joints that allow instrument levers to move laterally without causing wrist fatigue, every single ergonomic aspect is designed to optimise mobility.

The dentist's module is light and compact. Instruments are easy to use thanks to generous tubing excursion and a tubing recovery arm system minimises interference with the operating light. The handle incorporates a pushbutton to release the pneumatic lock on the module, helping the dentist position it perfectly after vertical adjustment.

All this ensures the dentist can always count on perfectly positioned instruments and, therefore, maximum working comfort.



Compact tray

The stainless steel tray is conveniently housed under the module. Positionadjustable, it can easily be used by either dentist or assistant. A compact tray of minimum size and bulk is available as an alternative.



Instrument excursion

Generous instrument tubing excursion ensures improved comfort in all operating positions around the patient.



Transthoracic position

The arm system has been designed to allow transthoracic rotation so the dentist's module can be positioned on the assistant's side.



SideFlex instrument levers

Removable and equipped with swing joints that also allow lateral movement, the levers alleviate wrist fatigue and are self-balancing. SideFlex instrument levers are available as an

AN ORDERLY WORK SPACE



Everything within reach to optimise organisation and productivity

Thanks to the broader arm system on the hanging-tube instruments model, the dentist's module can reach all the usual working positions around the patient. Made of light cast aluminium, the fixed arm can pass over the patient whatever the situation. A broad handle running the full width of the module makes module movement even easier. The handle can be removed for cleaning purposes. When placed in their holders, instruments are supported securely and arranged diagonally to optimise ergonomics. The syringe is positioned vertically at the start of the handpiece holder, ensuring a good grip. A large easy-to-read LCD touch screen is positioned above the

instruments. Dentists can control patient chair movements, hygiene systems and other integrated devices. Where the i-XS4 micromotor is present the instrument control panel also provides a full complement of endodontic functions (including control of the apex locator, with data display during canal treatment).



Tray holder
Large transthoracic tray holder
(optional) mounted on the
International module.



Sixth instrument
The international dentist's module can
be completed with an optional sixth
instrument: camera or T-LED light.







Assistant's module

The module mounted on the double-articulated arm is height-adjustable and comes with three holders as standard. An optional 5-holder version is also available: this can be used to hold other instruments such as a video camera or T-LED light. The keypad incorporates controls that move the patient chair, switch on the operating light and operate the water-to-cup and cuspidor bowl rinse functions.

FOR EVERY OPERATING STYLE

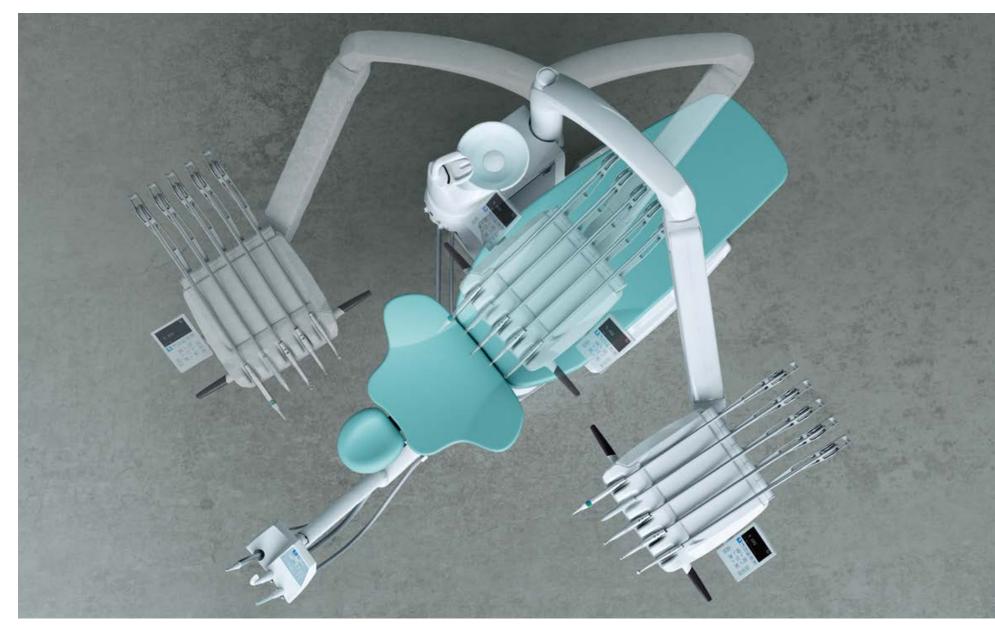



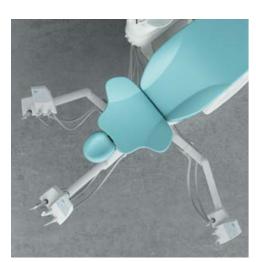

Ergonomic solutions that make the most of versatility

Equipped with a dentist's module that can be positioned anywhere around the treatment area, the Hybrid model features an assistant's module that can be shifted right or left, always remaining independent from dentist module rotation. Exceptional operating freedom whether working alone or with an assistant.


To complete module conversion where the unit is configured for left-handed dentists, the S220TR Hybrid has a repositionable instrument control panel: just detach it, rotate the quick coupling 180° and then reposition the panel on the opposite side of the module.

The S220TR Hybrid is suitable for all surgeries, even where spaces are narrow, and is available in the Continental version only.

Unmatched flexibility and high-strength materials

Compact and made of materials that make the dental unit structurally stable, the S220TR Side Delivery retains the geometry of the International module. Fitted on the cast aluminium column, the dentist's module can be height-adjusted via the control panel. Access to the patient chair is smooth and comfortable and, prior to treatment, the instruments remain outside the patient's field of view.

ENDODONTICS AT YOUR FINGERTIPS

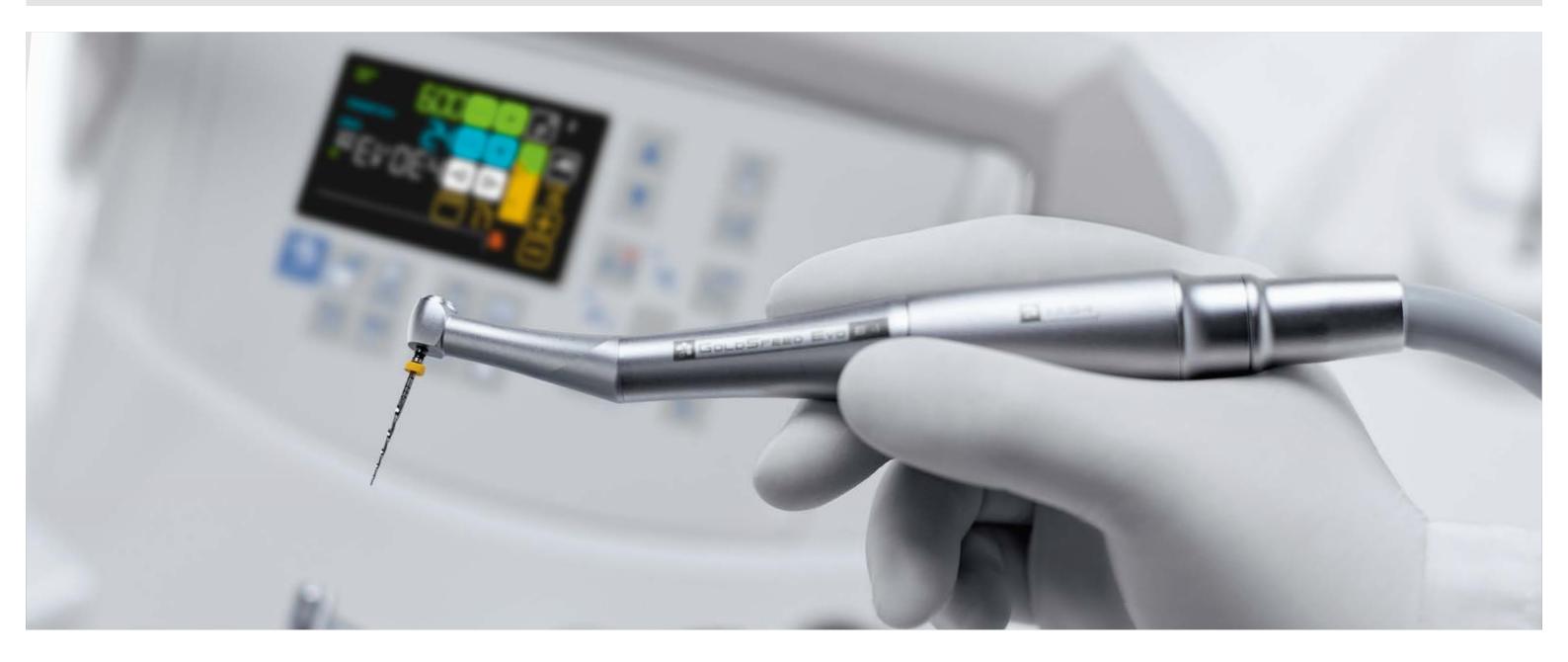

Instruments and software to perform canal treatment in complete safety

With the optional i-XS4 motor, used with the dedicated contra angle, the S220TR has everything needed for endodontic treatment thanks to fine adjustment of endocanular burr speed and torque.

An integrated apex locator is also available. During treatment, essential data is displayed in a clear user-friendly manner.



SC Scalers

Ultrasound scalers, compatible with the best commercially available tips, for supragingival prophylaxis and periodontal tasks. Also suitable for the preparation of small cavities and usable in endodontic root canal cleaning techniques. Available with or without LED lighting.



Micromotor i-XS4

Supplied as an alternative to the i-XR3 LED micromotor, the i-XS4 is suitable for conservative, prosthetic and endodontic dentistry. Speed ranges from 100 to 40,000 rpm and torque is adjustable. Powerful, light and versatile, it features autoforward and auto-reverse functions and is fully autoclavable.




Apex locatorWith the dedicated contra angle, the apex locator allows the apex distance to be displayed during the root canal instrumentation phase. The nearing of the apex is also indicated by a warning buzzer; to enhance safety, micromotor rotation is automatically disabled when it is reached.



EVO E4 contra angle

Featuring a 4:1 ratio specific for endodontics, the EVO E4 contra angle can be autoclaved and heatdisinfected. The miniaturised head makes access to the treatment zone easier. Recommended where the integrated apex locator is used.

MADE-TO-MEASURE COMFORT

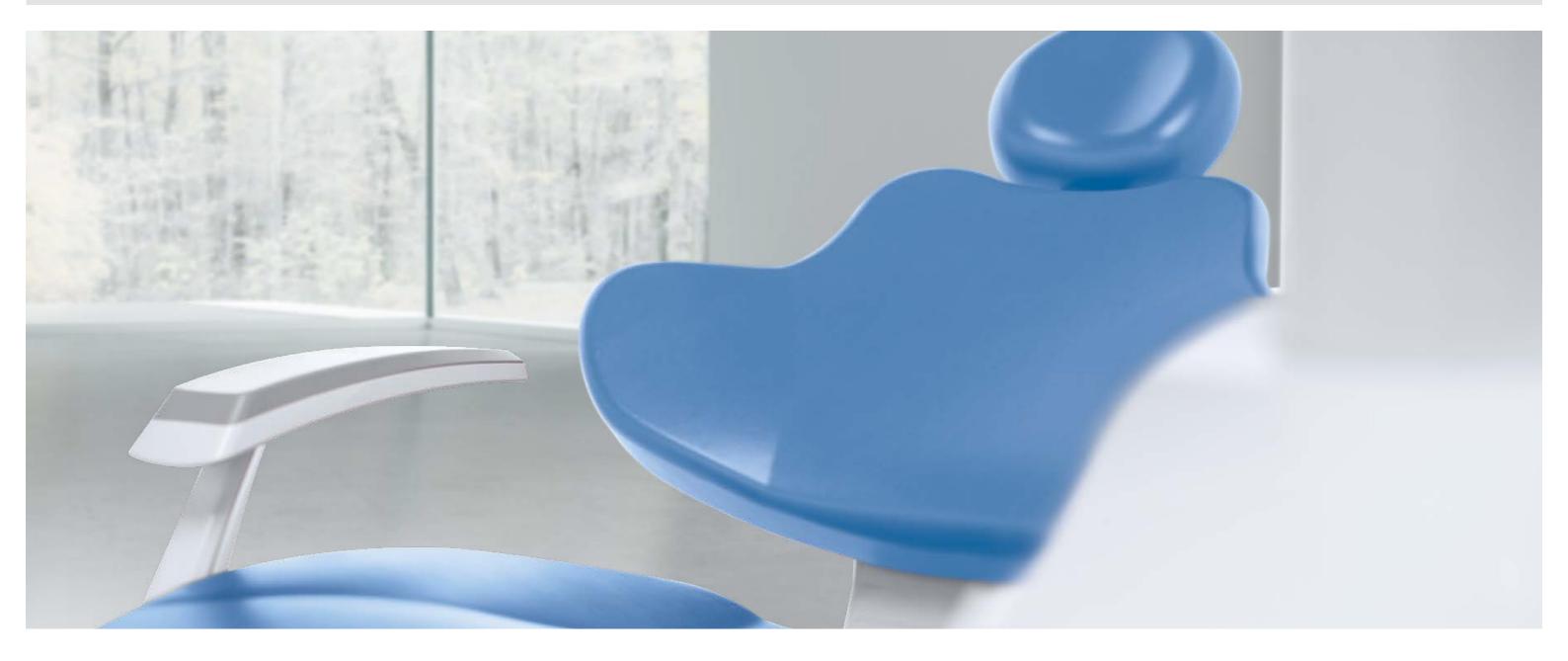

Designed for outstanding comfort

Alongside ergonomics, the comfort factor is of prime importance. Indeed, a patient at ease enhances the working effectiveness of the medical staff. Lined with hard-wearing upholstery in a wide choice of colours, the patient chair is deliberately shaped to ease the approach to the work area. Made of tough materials and featuring thoroughly tested mechanics, the patient chair offers excellent support and has a lift capacity of 190 kg.

The optional Nordic backrest ensures

lasting comfort during protracted treatment sessions and is shaped to aid dentists who work in 'indirect vision' mode.



Exceptional comfort and proper anatomical support thanks to special

upholstery options.



Universal headrest Supplied as standard, the headrest rotates on two axes and can be positioned for adults and children alike.



Atlaxis headrest


The headrest that follows the patient's anatomy. Ultra-simple move closer to the patient during fingertip release of the pneumatic treatment. Suitable for adults, lock allows orbital movement to children and any treatment type. achieve perfect positioning, vertical adjustment included.



Standard backrest

Shaped so that dental staff can

VISION AND CONTROL

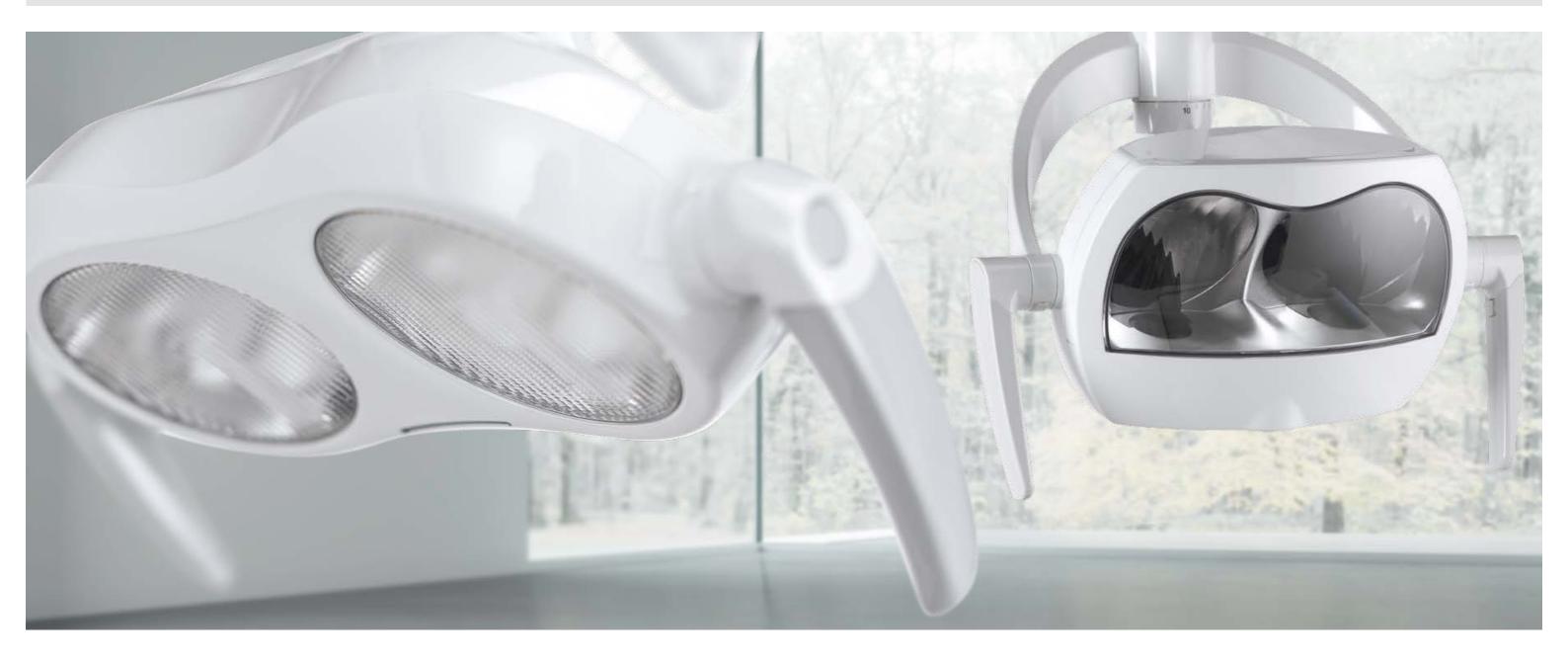

From lighting to foot controls, personal preferences matter

The optional Venus LED MCT is a major step forwards in the development of clinical lighting. The colour rendering index (CRI) exceeds 95, a factor that allows ultra-high colour fidelity. It's also possible to use 3 colour temperatures so users can select just the right lighting for their specific clinical disciplines. 4300K warm light for surgical treatment; 5000K neutral light for conservative work; 5500K cool light for realistic colour capture. Light intensity is adjustable up to 50,000 Lux. A broad spotlight and excellent depth of field (from 55 to 85 cm) limit the number of times the dentist has to reposition the lamp during treatment. What's more, "composite"

mode can be activated quickly to prevent premature composite polymerisation. Supplied as standard, the Venus Plus – L LED, with light intensity from 3,000 to 50,000 Lux and a colour temperature of 5000K, features three-axes rotation and incorporates the no-touch sensor for adjustment via potentiometer.



Multifunction foot control

Cutting-edge device with ergonomic lateral-slide control that lets dentists activate Chip Air and Chip Water functions on the instrument in use, invert micromotor rotation (when extracted), control patient chair movement and recall saved positions.



Pressure-operated foot control

Offering the same performance as the multifunction foot control, this unit offers different ergonomics for provides the same functions as the those who prefer a pressure-operated approach. Wireless versions of both this unit and the multi-function foot control are available: both have longlife batteries



Power Pedal

This third foot control model offers different design and ergonomics. It more classic models. The wireless Power Pedal version is not available.

HD MULTIMEDIA SYSTEM



Sharing information reassures patients and aids cooperation

Coupled with a latest-generation HD video camera, the Full HD 16: 9 tilt-adjustable 1920 x 1080 pixel flat screen medical monitor can be lead-connected to a PC. LED light sources ensure high levels of contrast and brightness while, thanks to the IPS panel which significantly broadens the screen viewing range, images can be viewed comfortably from any angle.



22" LED Monitor 93/42/EEC certified, the monitor is also available in the multitouch version; the screen can be positioned as desired.



C-U2 HD camera
The video camera

The video camera aids dentist-patient communication and, thanks to slender handpiece design, can reach distal zones with ease. Excellent depth of field avoids any need for manual focusing. Speeds up diagnostics and generates material with which to document treatment and complete clinical reports.


HD



Macro Cap

An extremely useful accessory, Macro Cap can magnify up to 100x in ultrahigh resolution. Three additional, high purity glass lenses optimise lighting of details located close to the optical unit.

REAL-TIME X-RAY DIAGNOSIS

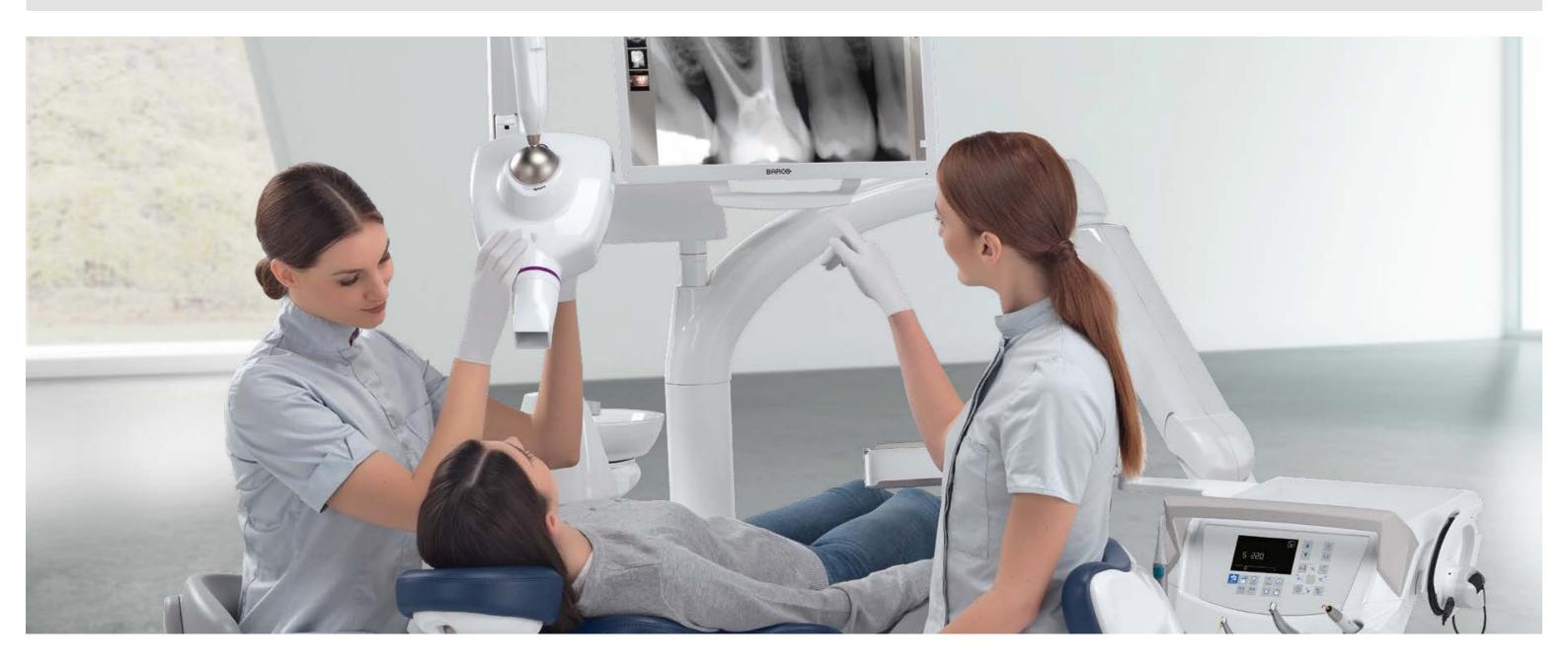

Unified in a single area, the system features cuttingedge X-ray devices Integrating a complete system boosts diagnostic potential while shortening treatment times. Featuring advanced technology, the system includes the X-ray unit with digital wireless handheld controller, a LED monitor with optional touch-screen function and a latest-generation digital sensor. X-ray images can be displayed in a matter of seconds on the 22" monitor, letting dentists share information with patients immediately before deciding on a treatment plan.



RXDC - HyperSphere technology

Incorporated on the dental unit, the X-ray unit and handheld control unit require no additional installation points. Thanks to a combination of 30 cm collimation and a 0.4 mm focal spot, parallelism is excellent and images are pin sharp. Locked and unlocked by touch-sensitive pads, the tube head rotates freely around the spherical joint.





Zen-X X-ray sensor

Extractable sensor with USB lead located on the dentist's module, can capture high-resolution images with minimal radiation dosage. Able to be sanitised, the sensor comes in two different sizes and is IP67 certified against water and dust infiltration.

COMPLETE PROTECTION



Meticulous design and certified devices ensure protection you can count on, always

Free to configure the dental unit with a broad selection of active devices as desired, dentists can also count on a design concept with an array of features that safeguard patient and surgery staff alike.



Activated via the control panel, the device ensures tubing is kept clean after the machine has remained idle for a time. Included in the BIOSTER or available separately.



Cuspidor bowl with optical sensor

The glass cuspidor bowl rotates and can be removed for cleaning purposes. Incorporated in the cup and bowl fill unit, the optical sensor automates cup filling.



A.C.V.S.

Automatic system that carries out suction system sanitisation cycles between one patient and another. Practical and fast.



The S220TR can be equipped with the automatic BIOSTER system. This device performs fast, intensive disinfection cycles on spray water circuits. Each phase of the cycle is automatically controlled by the dental unit. Dentists can personalise settings as desired. The cycle stages are:

• Introduction of compressed air into water pipes Withdrawal of disinfectant from tank and introduction into spray water circuits

• Presence of liquid in ducts for the required time • Liquid draining with compressed air • Flushing with (mains or distilled) water Tests performed by the Sapienza University of Rome and the Department of Public Health and Paediatrics of the ${\bf University\ of\ Turin\ }$ have shown that combining execution of a daily **BIOSTER** intensive disinfection cycle with use of the W.H.E. system results in the complete absence of bacterial load in the cooling liquids delivered by the

instruments.

Continuous Peroxy Ag+ (H₂O₂) based DIN disinfection system (compliant with **DVGW** EU standards and DVGW certified) and air gap separation device between mains water and dental unit circuits. Acts against all waterborne contaminants, including Legionella.



SANASPRAY

Tank with distilled water as an alternative to mains water



Removable International handle

The handle can be removed for in-depth sanitisation: removal also makes it easier to clean the dentist's module and instrument control



Removable Continental handle

New dentist's module handle can easily be removed for in-depth sanitisation.

SIDE DELIVERY

S220TR | VERSATILITY AT ITS VERY BEST

TECHNICAL CHARACTERISTICS

S220TR	CONTINENTAL	INTERNATIONAL	SIDE DELIVERY	CART	HYBRID
YGIENE SYSTEMS					
IOSTER	*	*	*	*	*
C.V.S.	*	*	*	*	*
ANASPRAY	*	*	*	*	*
H.E.	*	*	*	*	*
F.C.	*	*	*	*	*
ITIST'S MODULE					
) touch screen	•	•	•	•	•
23 micromotor with fibre optics (100 – 40,000 rpm, 3.3 Ncm)	•	•	•	•	•
4 micromotor with f.o. (100 – 40,000 rpm, 5.3 Ncm)	*	*	*	*	*
ler (with and without f.o.)	*	*	*	*	*
y syringe	•	•	•	•	•
instrument	*	*	*		
rated X-ray sensor	*	*	*	*	*
rated X-ray sensor	*	*	*	*	*
BODY	*	τ	*		Ť
cuspidor bowl	•	•	•	•	
nic cuspidor bowl					
dor bowl drive system	*	*	*	*	*
l sensor	*	*	*	*	*
ller water heater	*	*	*	*	*
	*	*	*	*	*
heating	*	*	*	*	*
ndent cannulae selection	*	*	*	*	*
stop device on chair base	•	•	•	•	•
nedia applications wiring	*	*	*	*	*
nt's tray holder	*	*	*	*	*
nt's module with 3 holders	•	•	•	•	•
ant's module with 5 holders	*	*	*	*	*
Plus LED lamp	•	•	•	•	•
LED MCT lamp	*	*	*	*	*
T CHAIR					
headrest	*	*	*	*	*
able right armrest	*	*	*	*	*
rmrest	•	•	•	•	•
s cushion	*	*	*	*	*
-function foot-control	•	•	•	•	•
ure-operated foot control	*	*	*	*	*
r Pedal foot control	*	*	*	*	*
ess multi-function foot-control	*	*	*	*	*
less pressure-operated foot control	*	*	*	*	*
	192 184 Memory Memory	Memory Me	195 182 emory Memory	Memory Me	197 187 emory Memory
oam Foam Foam Foam	Foam Foam	Foam F	oam Foam	Foam F	oam Foam
102 106 113 136	132 134	135	115 103	123	101 137
antic Mediterranean Pacific Indian B	132 134 lueberry Japanes		cottish Nevada	Polynesian C	101 137 aribbean Bright

INTERNATIONAL CONTINENTAL HYBRID CART